

*A*lotex[®]

 **Purgo** Dental
Biologics
Solution



“Feel the Clinical Freedom on Science and Safety”

Founded in 1999, Purgo Biologics strives to become one of the leading global companies in oral health care with its focus on safe biomaterials for soft tissue and bone regeneration.

Based on the specialized experiences accumulated by our outstanding research personnel, Purgo Research and Development Center based in Seoul is thriving to become the best in the world, specifically in the expertise of oral biomaterials for soft tissue and bone regeneration. All members in Research and Development Center are pursuing the optimized technical developments with various clinical studies, cooperative research with the governments, clinicians and educational institutions.

The solutions manufactured by Purgo are gaining fame throughout the world and Purgo’s solutions are widely accepted by global dentists from more than 30 countries.

Our production site is complying with the most international quality standards and regularly inspected by international agencies. Each production stage of our biologics solutions are controlled from the selection of the raw material to the final product.

Availability of words may vary from country to country.

Minimal memory, No tangle, and Superior handling



Biotex® Non-Resorbable PTFE Suture is comprised of a single-arm, non-resorbable monofilament suture with a stainless-steel surgical needle connected to the suture. The suture is uncoated, undyed and sterile for single use only, composed of 100% PTFE.

- ✔ SOFT HANDLING
- ✔ BIOLOGICALLY INERT
- ✔ NO TANGLE
- ✔ EASY KNOTTING





Indications

- Bone grafting procedures
- Periodontal surgery
- Guided tissue regeneration
- Ridge augmentation
- Implant surgery
- Soft tissue grafts

Specifications



Item N°	USP Size	Length (cm)	Needle Length (mm)	Circle	Point Type
BT3019	3-0	45	19	3/8	▽
BT3016	3-0	45	16	3/8	▽
BT4016	4-0	45	16	3/8	▽
BT4019	4-0	45	19	3/8	▽
BT5016	5-0	45	16	3/8	▽
BT4013	4-0	45	13	3/8	▽
BT5013	5-0	45	13	3/8	▽
BTP4013	4-0	45	13	1/2	⊙

Biotex™ Main Features & Benefits [10]

Suture

1 High pliability (PTFE)

- Tying and bending more at ease with less unintended loosening.


2 No room for little plaque.

- It dispels the possibility of any bacterial infection as well as the plaque formation and any other factors that prevent healing process.

Needle


1 Slim reverse cutting needle tip

- Precision slim cut triangular needle for small penetration area and smooth suturing.
- Minimize damage to surrounding soft tissue.

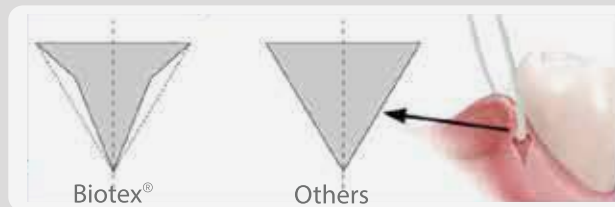


2 Strong Attachment

- Advanced technology for strong needle attachment .
- Smooth and firm connection between needle and thread .
- Rapid healing process due to the reduced bleeding from needle insertion .

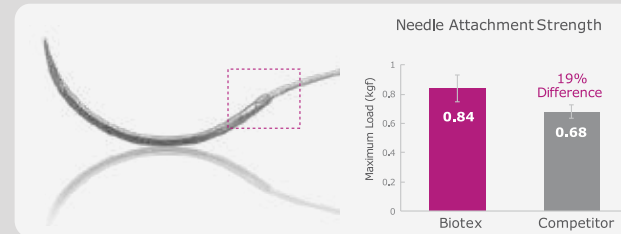

3 Strong Needle

- 33% higher strengths are required to bend needle in same degree compared to other product.
- High rigidity of the needle resists to bent stress during suturing.

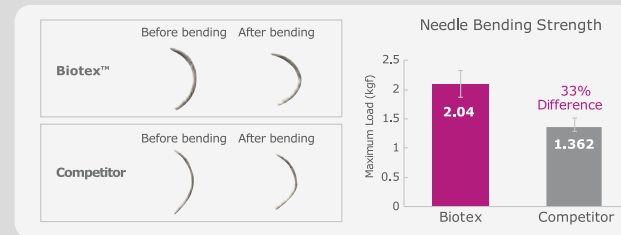

Needle Holding Clip

Designed to hold the needle in place, also allows for secure and easy release of the suture needle from its package.

Sturdy & Flexible Transparent Cover

Protect and give clear visibility of suture and needle.
Soft and sturdy cover effectively protect the suture.



Tab

Allows surgeon to easily grasp and remove the suture needle from its needle holder clip.

'Race Track' Shape

Designed to prevent suture from entangling, and allows easy release of the suture.



Benefits

- Soft and comfortable for patients
- Soft texture for patient comfort
- Reliable closure period
- Superior handling: provides flexibility in the positioning of a square knot. Easy to tie - Easy to remove
- Nonwicking: Elimination of bacterial wicking usually associated to monofilament
- Maintains tensile strength
- PFOA free

THE Graft™ 

THE Graft™  Collagen

OpenTex® 

OpenTex®-TR 

Biotex® 

BioCover™ 

 **Purgo** Dental
Biologics
Solution
Purgo Biologics Inc.

812, 27, Dunchon-daero 457beon-gil,
Jungwon-gu, Seongnam-si,
Gyeonggi-do, 13219, Korea
Tel: +82 2 548 1875
Made in Korea
www.purgo-biologics.com

EU Importer

Purgo Biologics Europe

1 Square Felix Bloch,
Pôle Activ' Ocean,
85300 Challans, France
Tel: +32 (0)2 28 10 61 02
E-mail: europel@purgobiologics.com
www.purgo-europe.com

EU Authorized Representatives

OBELIS S.A

Bd. Général Wahis, 53
1030 Brussels, Belgium
Tel: +32 2 732 59 54
Fax: +32 2 732 60 03
E-mail: mail@bbefe.net
www.obelis.net

THE Graft™, OpenTex®,
OpenTex®-TR, Biotex®,
THE Graft™ Collagen.